


MICROCNX™ SERIES DIMENSIONS AND WEIGHTS

POLYCARBONATE with teal cover For autoclave or gamma irradiation applications.

TERMINATION	METRIC EQ.	PART NO.	A	B	C	PART WEIGHT
1/8" HOSE BARB	3.2 mm	CNX17002HT	0.93 (23.6 mm)	2.25 (57.2 mm)	1.49 (37.8 mm)	0.27 oz (7.64 g)
3/32" HOSE BARB	2.4 mm	CNX17003HT	0.93 (23.6 mm)	2.16 (54.9 mm)	1.49 (37.8 mm)	0.27 oz (7.62 g)
1/16" HOSE BARB	1.6 mm	CNX17001HT	0.93 (23.6 mm)	2.05 (52.1 mm)	1.49 (37.8 mm)	0.27 oz (7.60 g)

PRODUCT DIMENSIONS

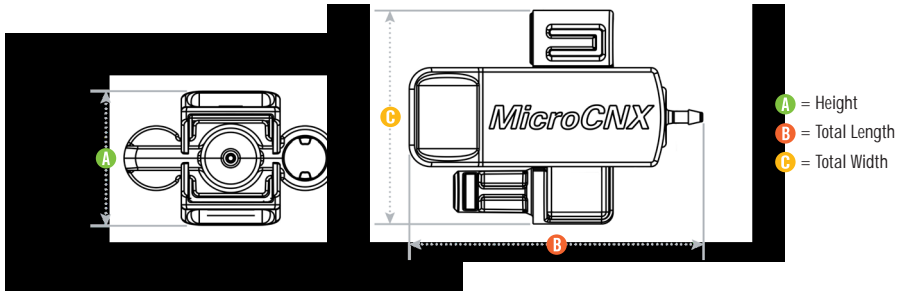

MICROCNX™ CONNECTORS ASSEMBLY PROCEDURE

PINCH

Remove the protective cover on each half.

CLICK

Join halves and click together.

PULL

Pull membrane strips directly away from the connector.